

CA 20:1, CA 20:1 L

CA 20:1 L KM

CA 20:1 L Micro-Series

CA 20:1 L KM Micro-Series



#### **ENG INSTRUCTIONS FOR USE**

Other languages available on https://dental.bienair.com/IFU



#### **Devices**













REF 1600632-001

REF 1600598-001

REF 1600692-001

REF 1600785-001

REF 1600786-001





REF 1500984-010



IRRIGATION LINE KM (10/pkg) REF 1501635-010



IRRIGATION LINE SYST KM 20:1 L (10/pkg) REF 1501621-010



MAINT SPRAYNET® (BOX of 6 CANS) REF 1600036-006



MAINT LUBRIFLUID® (BOX of 6 CANS) REF 1600064-006



MAINT AQUACARE (BOX of 6 CANS) REF 1600617-006

# Table of content

Cumphala	5.2 Changing the tool
. Symbols 4	5.2 Changing the tool
1.1 Description of symbols used 4	6. Maintenance and Servicing 14
2. Identification & Intended Use 5	6.1 Maintenance - General information 14
2.1 Identification 5	<b>6.1.1.</b> Precautions for maintenance 14
2.2 Intended Use 5	6.1.2. Suitable maintenance products 14
2.3 Intended patient population 5	6.2 Cleaning
2.4 Intended user 5	6.3 Disinfection
2.5 Use Environment	<b>6.3.1.</b> Manual cleaning and disinfection 16
2.6 Intended Medical Conditions 5	<b>6.3.2</b> .Automatic disinfection
2.7 Patient contra-indication and side effects . 5	6.4 Lubrication
2.8 In case of accident	6.5 Sterilization
8. User and Patient Safety: Warnings &	<b>6.5.1</b> .Procedure
Precautions for use	6.6 Packing and storage
Description 9	6.7 Servicing
4.1 Overview 9	7. Transport and Disposal
4.2 Technical Data	7.1 Transport
4.3 Classification	7.2 Disposal
4.4 Performances	8. General Information
4.5 Operating conditions	8.1 Terms of guarantee
5. Operation	8.2 References
5.1 Pictograms used	
i ictograffis useu	

## **ENG INSTRUCTIONS FOR USE**

## 1 Symbols

## 1.1 Description of symbols used

Symbol	Description Description	Symbol	Description
<b>C</b> € 0123	CE Marking with number of the notified body.	•••	Manufacturer.
$\triangle$	WARNING: hazard that could result in serious injury or damage to the device if the safety instructions are not correctly followed.	<u> i</u>	Consult instructions for use or consult electronic instructions for use.
<u> </u>	CAUTION: hazard that could result in light or moderate injury or damage to the device if the safety instructions are not correctly followed.		Data Matrix code for product information including UDI (Unique Device Identification).
X. X.	Temperature limitation.	% Trx	Humidity limitation.
,y	Atmospheric pressure limitation.	MD	Medical device.
<del>*</del>	Keep away from rain.	SN	Serial number.
	Wear rubber gloves.	REF	Catalogue number.
Rx Only	Warning: in accordance with federal law (USA), this device is only available for sale upon recommendation by an accredited practitioner.	EC REP	Authorized EC Representative in the European Community.
	General symbol for recovery/ recyclable.	135°C	Sterilization up to the specified temperature.
「」	Can be processed in an automated washer/disinfector for thermal disinfection.	-\\\\\\\	Lamp; lighting, illumination.

## 2 Identification & Intended Use

#### 2.1 Identification

Medical devices manufactured by Bien-Air Dental SA.

#### Type

Dental contra-angles handpiece (CA), push-button bur locking, with or without light, speed reducing ratio and external irrigation.

#### Description

Bien-Air Dental contra-angles are designed to transmit and apply the mechanical energy produced by an electric micromotor.

Contra-angle	Light		E-type connection (I	ISO 3964)
Ratio	With light	Without light	Standard	Micro-Series
•		•	•	
•	•		•	
•	•			•
•	•		•	
•	•			•

#### 2.2 Intended Use

Device intended for use in implantology.

#### 2.3 Intended patient population

The intended patient population for the contra-angles includes any person visiting a medical practitioner's office to receive treatment in line with the intended medical condition. There is no restriction concerning subject age, race, or culture. The intended user is responsible to select the adequate device for the patient according to the specific clinical application.

#### 2.4 Intended user

Device intended for professional use only. Used by dentists, dental professionals and oral surgeons.

#### 2.5 Use Environment

Professional healthcare facility environment.

#### 2.6 Intended Medical Conditions

Dental implantology is the treatment to replace one or more missing teeth.

#### 2.7 Patient contra-indication and side effects

No specific patient contra indication, side effects nor warnings exist for the contra-angle devices when the devices are used as intended.

# ENG

#### 2.8 In case of accident

If an accident occurs, the device must not be used until repairs have been completed by a qualified, authorized and trained technician in a repair center.

If any serious incident occurs in relation to the device, report it to a competent authority of your country, as well as the manufacturer through your regional distributor. Observe relevant national regulations for detailed procedures.

#### **⚠** WARNING

Any use other than that for which this device is intended is prohibited and may prove dangerous.

# 3 User and Patient Safety: Warnings & Precautions for use

#### **↑** WARNING

The device must be used by qualified dental professionals in compliance with the current legal provisions concerning occupational safety, health and accident prevention measures, and these instructions for use. In accordance with such requirements, the operators:

- Must only use devices that are in perfect working order; in the event of irregular functioning, coolant failure, excessive vibration, abnormal heating, unusual noise or other signs that may indicate malfunction of the device, the work must be stopped immediately; in this case, contact a repair center that is approved by Bien-Air Dental SA and have the service personnel carry out repair work.
- Must ensure that the device is used only for the purpose for which it is intended, must protect themselves, their patients and third parties from any danger.
- Any modification of the medical device is strictly forbidden.

#### To prevent any risk of infection, the warnings below must be observed:

#### **⚠** WARNING

- The device is supplied not sterile. To avoid any infection, respect the cleaning, disinfection, sterilization and overall maintenance procedure detailed in section 6. Sterilization before first use is mandatory.
- Medical personnel using or performing maintenance on medical devices that are contaminated
  or potentially contaminated must comply with universal precautions, in particular the wearing of
  personal protective equipment (gloves, goggles, etc.). Pointed and sharp instruments should be
  handled with great care.
- When disposing of the device, the user must return it sterilized to their dealer or contact an authorized body for the treatment and recovery of this type of equipment.
- While performing surgical treatment, the handpiece must not receive pressurised cooling air from the unit, to prevent contamination of the area being treated.

#### To prevent any risk of injury the warnings below must be observed:

#### **⚠** WARNING

- Never partially insert a bur to increase its active length.
- Do not touch the bur while the device is rotating.
- Each time a bur is inserted, check that the bur is fully inserted to the stop. Always check that the bur is locked by gently pushing and pulling the bur.
- Never use a bur if the tip is not compliant with specifications.
- Always ensure that the coolant supply is sufficient and adequate according to specifications.
- Soft tissues (tongue, cheeks, lips, etc.) must be protected by moving them away using a retractor or dental mirror to avoid the risk of burns if the push-button is inadvertently pressed while the instrument is running.

#### To prevent any risk of device overheating, the cautions below must be observed:

#### ⚠ CAUTION

- The device must not be started without a bur inserted into the chuck.
- To avoid overheating of the push-button, it should not be pressed inadvertently while the instrument is rotating.

#### To prevent any risk of device failure or malfunction the cautions below must be observed:

#### **A** CAUTION

- Before performing any clinical application, always test your device without any load to ensure it is in perfect order.
- Only use original Bien-Air Dental SA devices and accessories or those recommended by Bien-Air Dental SA.
- Respect the cleaning, sterilization and maintenance procedure detailed in section 6.
- These medical devices are intended to be used at a maximum height of 1.5 m, to avoid damage in the event of a fall.
- Never insert or remove a device while the micromotor is rotating.

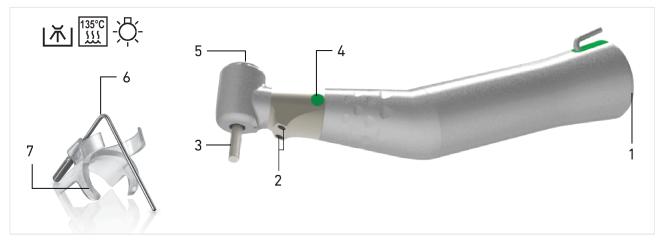


FIG. 1

## 4 Description

#### 4.1 Overview

#### FIG. 1

- (1) Micromotor connection
- (2) Light output (if applicable)
- (3) Tool (not supplied)
- (4) Transmission
- (5) Push-button with a bur-locking system
- (6) Kirshner/Meyer irrigation system
- (7) Plastic holder for the Kirshner/Meyer irrigation system

**Note**: The technical specifications, illustrations and dimensions contained in these instructions are given merely as an indication. They may not give rise to any claim.

The original language of these instructions for use is English.

For any further information, please contact Bien-Air Dental SA at the address given on the back cover.

#### 4.2 Technical Data

ii iooiiiioat Data	
Technical data	CA 20:1
Motor coupling compatibility	Coupling according to ISO 3964 - MS & MICRO-SERIES can be coupled to short and extra short Motor coupling - Other CA & PM can be coupled with all coupling type
Transmission ratio according to ISO 14457	Speed decreasing ratio 20:1 (green color)
Motor max speed	40'000rpm
Tool max speed	2'000 rpm
Irrigation type	External irrigation
Maximum Irrigation flow*	110 ml/min

<sup>\*</sup>Using Bien-Air Chiropro device at maximum irrigation level (5 drops).

#### ⚠ WARNING

- Always ensure that the coolant supply is sufficient and adequate according to specifications.
- The minimum irrigation quantity at the entrance of the CA must be set to at least 60mL/min (corresponding to 2 drops level if used in association with Chiropro devices).

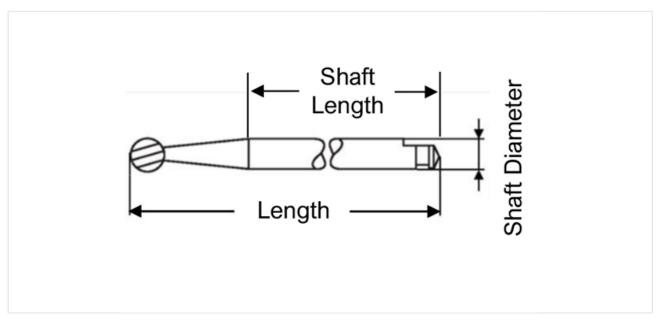


FIG. 2

#### Compatibility of the chuck mechanism

Shaft diameter 2.35 mm (0.09 in), type 1 as per ISO 1797, max. length 34 mm (0.87 in), code 4-6 as per ISO 6360-1 FIG. 2.

In implantology, longer cutting tools (Min. shaft length 12 mm in accordance with ISO 1797) can be used under the responsibility of the practitioner.

For implantology tools, minimum hardness is 50 HRC (520 HV) to prevent their deformation and jamming in the chuck mechanism.

Bur compatibility	
Shaft diameter ISO 1797	2.35 mm (Type 1)
Shaft length ISO 1797	≥ 12 mm
Bur max. working diameter	≤ 4 mm
Bur max. length ISO 6360-1	≤ 34 mm (code 4-6)
Implantology tool compatibility	
Shaft diameter ISO 1797	2.35 mm (Type 1)
Drill diameter	≤ 4.2 mm
Implant driver diameter	≤ 5.0 mm
Total length	≤ 41 mm

#### **MARNING**

- Never use a tool if the tip is not compliant with specifications.
- Follow the guidelines for use, according to bur manufacturer's instructions.

## 4.3 Classification

Class IIa in accordance with European Medical Regulation (EU) 2017/745.

## 4.4 Performances

Performances	CA 20:1
Speed transmission ratio	20:1

4.5 Operating conditions

Operating conditions		
x:x	Temperature range:	[+10°C; +35°C] [+50°F; +95°F]
X5_ X56	Relative humidity range:	[30%; 80%]
(+) • (+)	Air pressure range: Barometric pressure range:	[700 hPa; 1060 hPa] [490 mmHg; 795 mmHg]
<del>*</del>	Keep away from rain	

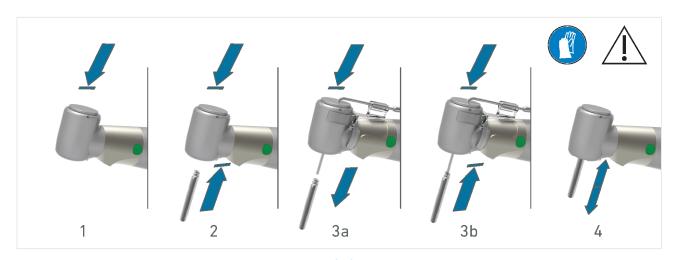


FIG. 3

## 5 Operation

#### 5.1 Pictograms used

Sym	Description	Sym	Description
$\bigcirc$	Movement in the direction indicated.	$Q_{\underline{\downarrow}}$	Movement to the stop in the direction indicated.
<b>†</b>	Back and forth		

### 5.2 Changing the tool

#### Clamping system operation

Push-button bur locking.

movement.

- 1. Press the push-button and simultaneously pull out the tool (FIG. 3 step 1).
- 2. Press the push-button, insert the new tool all the way to the stop and release the push-button (FIG. 3 step 2).
- 3. If the Kirschner Meyer irrigation system is used, it is recommended to leave the cooling tube in place to change the tool (FIG. 3 step 3a & 3b).

#### Clamping system inspection

Verify that the bur can be turned around smoothly and check that the bur remains in position when traction is applied (FIG. 3 step 4).

#### **⚠** WARNING

- Never insert or remove a bur while the device is rotating.
- Never partially insert a bur to increase its active length.
- Do not touch the bur while the device is rotating.
- Each time a bur is inserted, check that the bur is fully inserted to the stop. Always check that the bur is locked by gently pushing and pulling the bur.

#### **↑** CAUTION

- If the bur cannot be easily and fully inserted into the chuck, contact your usual supplier or Bien-Air Dental SA for repair.
- Test your device without any load to ensure the bur rotates stably and its dynamic eccentricity is acceptable for the planned clinical procedure.

## 6 Maintenance and Servicing

#### 6.1 Maintenance - General information

#### ♠ WARNING

- The device is supplied "nonsterile". Clean, dry, lubricate and sterilize the device prior to first use.
- Follow your national directives, standards and guidelines for cleaning and sterilization recommendations.
- Rest the device on a cleanable support to avoid risks of infection for yourself, the patient or third parties.

#### 6.1.1 Precautions for maintenance

- Within a maximum of 30 minutes after each treatment, clean and disinfect the instrument. Observing this procedure eliminates any blood, saliva or residues and prevents the transmission system from being blocked.
- Only use original Bien-Air Dental SA maintenance products and parts or those recommended by Bien-Air Dental SA. For suitable maintenance products refer to section Suitable maintenance products. Using other products or parts may cause faults during operation and/or void the guarantee.
- The device can be used for more than 1,000 sterilization cycles.
- The irrigation system can be used for 15 reprocessing (cleaning, disinfection and sterilization) cycles before disposal.

#### 6.1.2 Suitable maintenance products

#### Preliminary cleaning:

- Use tap water if the local tap water has pH within the range 6.5 8.5 and chloride content below 100 mg/l. If the local tap water does not meet these requirements, use demineralized (deionized) water instead.
- Aquacare.

#### Manual cleaning:

• Spraynet®.

#### Manual disinfection:

Alkaline detergent, or detergent-disinfectant (pH 8-11) recommended for cleaning-disinfection
of dental or surgical instruments. Disinfectant products composed either of
didecyldimethylammonium chloride, quaternary ammonium carbonate or neutral enzymatic
product. (e.g. Neodisher® mediclean) are also allowable.

#### Automatic cleaning-disinfection:

• Use an alkaline or enzymatic product recommended for cleaning in a washer-disinfector for dental or surgical instruments (pH 8-11).

#### Lubrification:

Lubrifluid®.





#### 6.2 Cleaning

- Do not use saline solution to keep the device moist until it can be cleaned.
- Clean using manual cleaning or automated washer/disinfector only (do not use ultrasonic cleaner).
- Carry out cleaning and sterilization without a bur in the chuck mechanism.

#### Preparation

- 1. Remove the irrigation line then disconnect the device from the electrical motor.
- 2. Remove the bur (FIG. 3 step 1).
- 3. If there is a large amount of debris, clean the exterior of the device with disinfectant wipes.

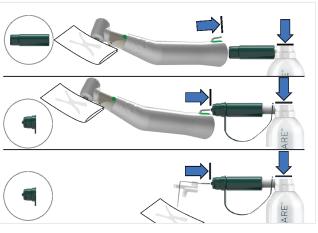


FIG. 5

#### Remove dirt / deposits

- 1. Clean the exterior and interior of the device under tap water at 15°C-38°C (59°F-100°F) provided that the local tap water has a pH within the range of 6.5 - 8.5 and a chloride content below 100 mg/l. If the local tap water does not meet these requirements, use demineralized (deionized) water instead. (FIG. 4)
- 2. After selecting the appropriate nozzle, perform preliminary cleaning of the device by using the product Aquacare. Spray the inside and the outside of the device and inside the irrigation tube (FIG. 5).

#### 

- Use the Aguacare maintenance pro duct as soon as possible, to remove NaCl deposits.
- If there is a large amount of debris, clean the exterior of the device with disinfectant wipes. Observe the instructions given by the manufacturer.



#### 6.3 Disinfection

#### 6.3.1 Manual cleaning and disinfection

- 1. Dip the device in a bath containing a cleaning and disinfectant product (e.g. an alkaline product like neodisher Mediclean). Follow the concentration and duration recommended by the fabricant of the disinfection product.
- 2. Brush the device with a smooth, flexible brush (e.g. soft-bristled toothbrush). DO NOT USE a wire brush.
- 3. Optional: perform additional cleaning and disinfection of the external surfaces with non-woven wipes impregnated with a disinfection product (e.g. didecyldimethylammonium chloride).
- 4. Rinse the exterior and interior of the device under tap water at 15°C-38°C (59°F-100°F) provided that the local tap water has a pH within the range of 6.5 8.5 and a chloride content below 100 mg/l. If the local tap water does not meet these requirements, use demineralized (deionized) water instead.
- 5. After selecting the appropriate nozzle, spray inside the device with Spraynet® (FIG. 6).
- 6. Dry the external surfaces with sterile nonwoven compresses (low linting textiles), preferably impregnated with Spraynet® or other blends of drying alcohols, like ethanol or isopropyl alcohol suitable for metals and polymers.

#### FIG. 6

#### 6.3.2 Automatic disinfection

**Note**: The automatic cleaning-disinfection can replace the previous steps 4 to 6.

#### Washer-disinfector

Carry out automatic cleaning- disinfection using an approved washer-disinfector which complies with ISO standard 15883-1.

#### Detergent and washing cycle

Use cleaning products (e.g alkaline de tergent or detergent-disinfectant pH

8-11 or neutral enzyme detergent pH 7-8) recommended for washer-disinfector.

#### Recommended specifications for the thermo-disinfection cycle.

Phase	Parameters
Pre-cleaning	<45°C (113°F); ≥ 2 minutes
Cleaning	45-55°C/113-131°F for enzymatic detergents and 45-65 °C/113-149°F for alkaline detergents ≥ 5 minutes
Neutralization	≥ 2 minutes
Rinsing	Tap water, ≤30°C (86°F), ≥ 2 minutes cold water
Thermal Disinfection	Demineralized water, 90°C — 95°C (194°F — 203°F), 5-10 minutes
Drying	18 — 22 minutes

Never rinse the devices to cool them.

If an automatic washer is used instead of the washer/thermo-disinfector, respect the previous program for the Pre-cleaning, Cleaning, Neutralization and Rinsing phases. If the local tap water has a pH outside the range of 6.5-8.5 or if it contains more than 100 mg/l chloride (Cl-ion), do not dry the device inside the automatic washer but dry it manually with low linting textiles.

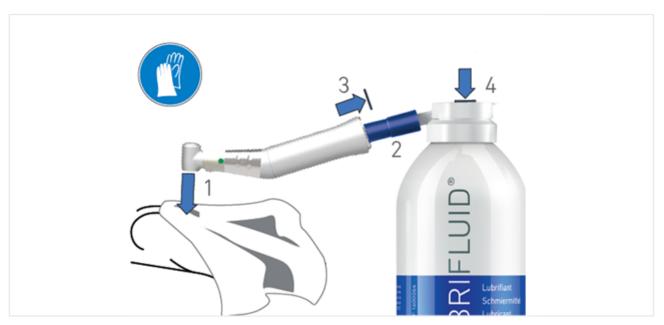


FIG. 7

#### 6.4 Lubrication

#### Verifying cleanliness

Visually inspect the device to ensure it is clean. Repeat the cleaning and disinfection procedure if necessary.

#### Lubrication

Lubricate before each sterilization or at least twice a day. Only the Lubrifluid® spray must be used. FIG. 7

- 1. Place the device in a sterile, non-woven cloth to collect the excess of lubricant.
- 2. Select the appropriate nozzle.
- 3. Insert the nozzle of the Lubrifluid® can in the rear of the device's handle.
- 4. Activate the spray for 1 second and clean the excess oil on the exterior with a sterile, non-woven compress.

#### 6.5 Sterilization

#### ⚠ CAUTION

The quality of the sterilization is highly dependent on how clean the instrument is. Only perfectly clean instruments may be sterilized.

- To improve the effectiveness of the sterilization, make sure the contra-angle is completely dry before and after the sterilization.
- Do not use a sterilization procedure other than the one described below.
- Only use dynamic air removal cycles: pre-vacuum or steam flush pressure pulse (SFPP) cycles.
- If the sterilization is required by national directives, use only dynamic sterilizers: do not use a steam sterilizer with a gravity displacement system.
- As with all instruments, following each sterilization cycle, including drying, remove the device to avoid an excessive exposure to heat which can result in corrosion.
- The bur must be disassembled from the device prior to sterilization.

#### 6.5.1 Procedure

- 1. Pack the device in a packaging approved for steam sterilization.
- 2. Sterilize using steam, following dynamic air removal cycle (ANSI/AAMI ST79, Section 2.19), i.e. air removal via forced evacuation (ISO 17665-1, ISO/TS 17665-2) at 135°C (275°F), for 3 minutes or 132°C (269.6°F) for 4 minutes. In jurisdictions where sterilization for prions is required, sterilize at 135°C (275°F) for 18 minutes.

#### The recommended parameters for the sterilization cycle are:

- The maximum temperature in the autoclave chamber does not exceed 137°C (278°F), i.e. the nominal temperature of the autoclave is set between 132°C 135.5°C (269.6°F 275.9°F) taking into account the uncertainty of the sterilizer as regards temperature.
- The maximum duration of the interval at the maximum temperature of 137°C (278 °F) is in accordance with national requirements for moist heat sterilization and does not exceed 30 minutes.
- The absolute pressure in the chamber of the sterilizer is comprised in the interval between 0.07 bar to 3.17 bar (1 PSIA to 46 PSIA/28 Hg to 31 PSIG).
- The rate of change of temperature does not exceed 15°C/min (27°F/min) for increasing temperature and -35°C/min (-63°F/min) for decreasing temperature.
- The rate of change of pressure does not exceed 0.45 bar/min (6.6 psia/min) for increasing pressure and -1.7 bar/min (-25 psia/min) for decreasing pressure.
- No chemical or physical reagents are added to the boiler feedwater.

#### 6.6 Packing and storage

Storage conditions		
***	Temperature range:	[0°C; +40°C] [+32°F; +104°F]
,	Relative humidity range:	[10%; 80%]
	Air pressure range: Barometric pressure range:	[650 hPa; 1060 hPa] [490 mmHg; 795 mmHg]
<del>*</del>	Keep away from rain	

The device must be stored inside the sterilization pouch in a dry and dust free environment. The temperature must not exceed 55°C (131°F). If the device will not be used for 7 days or more after the sterilization, extract the device from the sterilization pouch and store it in the original package. If the device is not stored in a sterilization pouch or if the pouch is no longer sterile, clean, lubricate and sterilize the device before using it.

#### 

If the medical device has been stored refrigerated, allow it to warm up to room temperature prior to removal from its packaging and use.

#### A WADNING

Comply with the expiration date of the sterilization pouch which depends on the storage conditions and type of packaging.

#### 6.7 Servicing

Bien-Air Dental SA recommends a regular service for the handpiece after 4,000 processing cycles or 5 years. Never dismantle the device. For all servicing or repair operations, you are advised to contact your usual supplier or Bien-Air Dental SA directly.

# ENG

## 7 Transport and Disposal

## 7.1 Transport

Transport conditions		
х. х.	Temperature range:	[-20°C; +50°C] [-4°F; +122°F]
X5	Relative humidity range:	[5%; 80%]
	Air pressure range: Barometric pressure range:	[650 hPa; 1060 hPa] [490 mmHg; 795 mmHg]
<del>*</del>	Keep away from rain	

### 7.2 Disposal



The disposal and/or recycling of materials must be performed in accordance with the local, national or international regulations.

All contra-angles and handpieces must be recycled. In order to avoid any risk of contamination, the user must return the device sterilized to their dealer or contact an authorized body for the treatment and recovery of this type of equipment.

## 8 General Information

#### 8.1 Terms of quarantee

Bien-Air Dental SA grants the user a guarantee covering all functional defect, or material or manufacturing faults.

The guarantee period is 12 months from the date of invoicing.

In the event of a justified claim, Bien-Air Dental SA or its authorized representative will repair or replace the product free of charge.

All other claims of any kind whatsoever, particularly claims for damages, are excluded.

Bien-Air Dental SA cannot be held liable for damage or injury and the consequences thereof, resulting from:

- · Excessive wear and tear,
- Infrequent or improper use,
- Failure to observe the servicing assembly or maintenance instructions,
- Damage caused by unusual chemical, electrical or electrolytic influences.

#### ⚠ CAUTION

The guarantee becomes null and void if damage and its consequences result from incorrect servicing or modification by third parties not authorized by Bien-Air Dental SA. Guarantee requests will only be taken into consideration if the product is accompanied by a copy of the invoice or delivery note. The following information must be clearly indicated: purchase date, product reference and serial number.

#### 8.2 References

U.Z INCI	el elices
REF	Legend
1600598-001	Contra-angle handpiece speed reducing ratio CA 20:1 L with light, push-button grip, external irrigation
1600632-001	Contra-angle handpiece speed reducing ratio CA 20:1 without light, push-button grip, external irrigation
1600692-001	Contra-angle handpiece speed reducing ratio CA 20:1 L Micro- Series with light, push-button grip, external irrigation
1600785-001	Contra-angle handpiece speed reducing ratio CA 20:1 L KM with light, push-button grip, external irrigation + Kirschner Meyer type detachable irrigation system. Including irrigation set REF 1501621-010
1600786-001	Contra-angle handpiece speed reducing ratio CA 20:1 L KM Micro- Series with light, push-button grip, external irrigation + Kirschner Meyer type detachable irrigation system. Including irrigation set REF 1501621-010
1501621-010	Kirschner/Meyer type detachable irrigation set for CA 20:1 L KM and CA 20:1 L KM Micro-Series including 10 rings and 10 tubes
1501635-010	Kirschner/Meyer pack of 10 disposable sterile lines
1500984-010	Pack of 10 disposable sterile lines
1600617-006	Aquacare, cleaning spray for physiological liquid 500 ml, box of 6 cans
1600036-006	Spraynet®, cleaning spray 500 ml, box of 6 cans
1600064-006	Lubrifluid®, lubricant 500 ml, box of 6 cans



#### Bien-Air Dental SA

Länggasse 60 Case postale 2500 Bienne 6 Switzerland Tel. +41 (0)32 344 64 64 Fax +41 (0)32 344 64 91 dental@bienair.com

Other addresses available at www.bienair.com

EC REP Bien-Air Europe Sàrl
19-21 rue du 8 mai 1945
94110 Arcueil
France