

MICROMOTORS

- MX-i
- MX-I LED 3RD GEN
- MX-i LED



ENG INSTRUCTIONS FOR USE.

Other languages available on
<https://dental.bienair.com/IFU>

CE Rx Only
0123 REF 2100245-0010/2024.02

Packaging content (REF)



MOT MX-i
REF 1600825-001



MOT MX-i LED 3RD GEN
REF 1601008-001




MOT MX-i LED
REF 1600755-001

Optional accessories



CABLE MX-i LED
REF 1601069-001
(FOR MX-i and MX-i LED)

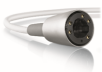

CABLE MX-i LED 3RD GEN
REF 1601009-001
(Only for MX-i and MX-i LED)



CABLE MX LED 3M
REF 1600881-001
(Only for MX-i LED 3RD GEN)



10X

O-RING 8.1x0.73
REF 1600967-010
(For all MX-i micromotors)



6X

MAINT SPRAYNET® (BOX OF 6 CANS)
REF 1600036-006








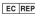









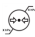



Table of contents

1 Symbols	4	5.5 Packing and storage	16
1.1 Description of symbols used	4	5.6 Servicing	16
2 Identification & Intended Use	5	6 Transport & disposal	17
2.1 Identification	5	6.1 Transport	17
2.2 Intended use	5	6.2 Disposal	17
2.3 Intended patient population	5	7 General information	18
2.4 Intended user	5	7.1 Terms of guarantee	18
2.5 Use Environment	5	7.2 Cable compatibility	18
2.6 Intended medical conditions	5	7.3 References	19
2.7 Patient contra-indications and side effects	5		
2.8 In case of accident	5		
User and Patient Safety: Warnings and Precautions for use:	6		
3.1 Installation	7		
4 Description	8		
4.1 Overview	8		
4.2 Assembly and preparation	9		
4.3 Technical data	10		
4.4 Classification	10		
4.5 Performances	11		
4.6 Operating conditions	11		
5 Maintenance and servicing	12		
5.1 Maintenance - General information	12		
5.2 Cleaning and disinfection	13		
5.3 Sterilization	15		
5.3.1 Procedure	15		
5.4 Lubrification	15		

ENG INSTRUCTIONS FOR USE

1 Symbols

1.1 Description of symbols used

Symbol	Description	Symbol	Description
	Manufacturer.		Catalogue number.
	CE Marking with number of the notified body.		Serial number.
	WARNING: hazard that could result in serious injury or damage to the device if the safety instructions are not correctly followed.		Medical device.
	CAUTION: hazard that could result in light or moderate injury or damage to the device if the safety instructions are not correctly followed.		Authorized EC Representative in the European Community.
	Wear rubber gloves.		Thermo washer disinfectable.
	Electrical security. Applied part type B.		Sterilizable in a steam sterilizer (autoclave) at the specified temperature.
	Warning: in accordance with federal law (USA), this device is only available for sale upon recommendation by an accredited practitioner.		Consult instructions for use or consult electronic instructions for use.
	Data Matrix code for product information including UDI (Unique Device Identification).		Temperature limit.
	Humidity limitation.		Atmospheric pressure limitation.
	Keep away from rain.		General symbol for recovery/recyclable.
	Recyclable electrical and electronic material.		

2 Identification & Intended Use

2.1 Identification

Medical devices manufactured by Bien-Air Dental SA.

Type:

MX-i LED micromotor

Brushless, self-ventilated, with light, fully sterilisable motors, machine washable

MX-i micromotor

Brushless, self-ventilated, without light, fully sterilisable motors, machine washable

MX-I LED 3RD GEN

Brushless, non-ventilated, with light, fully sterilizable motors, machine washable

These motors are powered through connection to dental unit / chair or console electronics. No spray features.

Description:

Bien-Air Dental micromotors are designed to transform electricity into mechanic rotation to drive dental straight hand-pieces and contra-angles.

2.2 Intended use

MX-I LED 3RD GEN

The device is intended for use in implantology.

MX-i and MX-i LED micromotor

Product intended for use in Implantology, Periodontology and Oral surgery.

2.3 Intended patient population

The intended patient population for the device includes any person visiting a dental practitioners' office to receive treatment in line with the intended medical condition. There is no restriction concerning subject age, race, or culture. The intended user is responsible to select the adequate device for the patient according to the specific clinical application.

2.4 Intended user

Product intended for professional use only. Used by dentists, dental professionals and oral surgeons.

2.5 Use Environment

Professional healthcare facility environment.

2.6 Intended medical conditions

- Dental implantology is the treatment to replace one or more missing teeth.
- Oral surgery treatments include impacted teeth extraction, wisdom teeth extraction, non-salvageable decayed teeth extraction, Guided and not-guided bone regeneration, apicoectomy, osteotomy, sequestrectomy and hemisection.
- The main periodontology treatments include gingivitis, periodontitis and surgical procedure therapy called flap surgery.

2.7 Patient contra-indications and side effects

No specific patient contra indication, side effects nor warnings exist for the device when it is used as intended.

2.8 In case of accident

If an accident occurs, the device must not be used until repairs have been completed by a qualified, authorized and trained technician in a repair center.

If any serious incident occurs in relation to the device, report it to a competent authority of your country, as well as the manufacturer through your regional distributor. Observe relevant national regulations for detailed procedures.



WARNING

Any use other than that for which this device is intended is prohibited and may prove dangerous.

3 User and Patient Safety: Warnings and Precautions for use:

WARNING

The device must be used by qualified dental professionals in compliance with the current legal provisions concerning occupational safety, health and accident prevention measures, and these instructions for use. In accordance with such requirements, the operators:

- Must only use devices that are in perfect working order; in the event of irregular functioning, coolant failure, excessive vibration, abnormal heating, unusual noise or other signs that may indicate malfunction of the device, the work must be stopped immediately; in this case, contact a repair center that is approved by Bien-Air Dental SA and have the service personnel carry out repair work.
- Must ensure that the device is used only for the purpose for which it is intended, must protect themselves, their patients and third parties from any danger.
- Any modification of the medical device is strictly forbidden.

Electrical safety:

WARNING

Electrical safety can only be claimed when the device is used with Bien-Air Dental compatible consoles, drives motors and cables.

- Always refer to the dental unit or console instruction for use to confirm compatibility with the device and electrical safety compliance.
- When used with Bien-Air dental medical devices, the device is compliant with standard IEC 60601-1 for which the following requirements apply:
 - Protection against electrical shock
 - Ingress of liquids
 - Protection against excessive temperatures and other safety hazards.

Electromagnetic compatibility:

WARNING

Electromagnetic compatibility can only be claimed when the device is used with Bien-Air Dental compatible consoles, drive motors and cables.

- Always refer to the dental unit or console instruction for use to confirm compatibility with the device and electromagnetic compliance.
- Since compliance with the international standard IEC 60601-1-2 does not guarantee immunity against 5G worldwide (due to the different frequency bands used locally), avoid the presence of devices equipped with 5G broadband cellular networks in the clinical environment or ensure that the network functionality of these devices is disabled during the clinical procedure.
- Magnetic disturbance can occur from other electromedical devices, refer to the dental unit or console instruction for use for EMC specifications.

To prevent any risk of explosion, the warnings below must be observed:

WARNING

According to IEC 60601-1:2005 +A1 2012 / AnnexG, electrified devices (motors, control units, couplers and attachments), can be safely used in a medical environment in which potentially explosive or flammable mixtures of anaesthetic substances are delivered to the patient only if:

- The distance between the motor and the anaesthetic breathing circuit exceeds 25 cm.
- The motor is not used simultaneously to the administration of the anaesthetic substances to the patient.

To prevent any risk of motor overheating, the warning below must be observed:

 **CAUTION**

The MX-I LED 3RD GEN is not self-ventilated. Do not use for applications outside of the intended use.

To prevent any risk of infection, the warnings below must be observed:

 **WARNING**

- The device is supplied not sterile. To avoid any infection, respect the cleaning, disinfection, sterilization and overall maintenance procedure detailed in section 5.
- Medical personnel using or performing maintenance on medical devices that are contaminated or potentially contaminated must comply with universal precautions, in particular the wearing of personal protective equipment (gloves, goggles, etc.). Pointed and sharp instruments should be handled with great care.
- Rest the device on a cleanable support.

To prevent any risk of injury and/or material damage the cautions below must be observed:

 **CAUTION**

- In the event of excessive vibrations, abnormal heating, unusual noise or other signs suggesting that the device is malfunctioning, work must be suspended immediately. In this case, contact a repair centre approved by Bien-Air Dental SA.
- Never connect an instrument on a running micromotor.
- Do not spray any lubricant or cleaning solution into the motor.
- Never rinse the device to cool it down.

3.1 Installation

 **WARNING**

The micromotor can only be used with Bien-Air Dental or compatible electronic devices. To conform with IEC 60601-1 safety requirement, only medical device power supply with two MOPP should be used.



FIG. 1

4 Description

4.1 Overview

FIG. 1

- (1) Motor nose
- (2) Motor body
- (3) Cable/motor connection

Nota : The technical specifications, illustrations and dimensions contained in these instructions are given merely as an indication. They may not give rise to any claim.

The original language of those instructions for use is English.

For any further information, please contact Bien-Air Dental SA at the address given on the back cover.

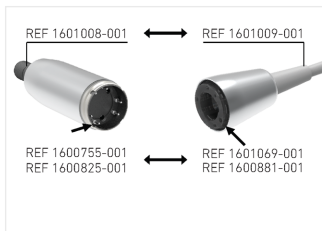


FIG. 2

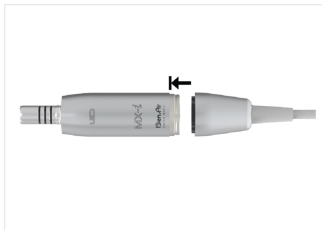


FIG. 3



FIG. 4



FIG. 5

4.2 Assembly and preparation

1. The MOT MX-i REF 1600825-001 & MOT MX-i LED REF 1600755-001 motors are installed using the REF 1601069-001 or REF 1600881-001 power cable. **FIG. 2.**
2. The MOT MX-i LED 3RD GEN REF 1601008-001 motor is installed with the REF 1601009-001 mains cable as shown in **FIG. 2.**
3. Check that the rear of the motor and the cable connector are clean and dry.
4. Position the motor and its proprietary cable as shown in **FIG. 2.**
5. The locator must be in front of the hold. Rotate it to find the exact position and push. **FIG. 3.**
6. Tighten the nut cable sleeve to the stop. **FIG. 4.**

Changing the seals manually (no tool required) FIG. 5

CAUTION

- Do not lubricate the O-ring.
- Use proprietary O-ring only.
- Check that the O-rings are neither broken nor scratched after mounting them.

4.3 Technical data

Technical data	
Coupling	Nose in accordance with ISO 3964*.
Operating times	No limitations for the user. Operating times are electronically imposed by Bien-Air control boards, as a function of the applied torque.
Rotation speed range	100-40'000 rpm
Direction of rotation	Clockwise and anticlockwise.
Luminous intensity	LED, 14 klux**
Maximum torque	800 mNm***

*Compatible with handpieces with dimensions "short", "middle" and "long".

**Measured in combination with REF1600632-001 CA 20:1

***Torque available at the tool when the motor is used in implantology with the CA 20:1.

4.4 Classification

Class IIa in accordance with European Medical Regulation (EU) 2017/745.

Class II type B device in accordance with IEC 60601-1 standard.

4.5 Performances

Performances

Give speed and torque as preset	Electronically controlled
Speed value accuracy	$\pm 5\%$

4.6 Operating conditions

Operating conditions



Temperature range:

$+10^{\circ}\text{C} - +35^{\circ}\text{C}$ ($+50^{\circ}\text{F} - +95^{\circ}\text{F}$)



Relative humidity range:

30% — 80%



Air pressure range:

700 hPa — 1060 hPa



FIG. 6

5 Maintenance and servicing

5.1 Maintenance - General information

Clean, disinfect, dry and sterilize the device prior to first use.

Within a maximum of 30 minutes after each treatment, clean the motor. Observing this procedure eliminates any blood or saliva residues.

 **CAUTION**

Follow your national directives, standards and guidelines for cleaning and sterilization recommendations.

 **CAUTION**

Do not spray any lubricant or cleaning solution into the motor [FIG. 6](#).

Suitable maintenance products:

Only use original Bien-Air Dental SA maintenance products mentioned below and parts or those recommended by Bien-Air Dental SA. Using other products or parts may cause faults during operation and/or void the warranty.

- Spraynet®
- Alkaline detergent or detergent disinfectant (pH 8-11) recommended for cleaning-disinfection of dental or surgical instruments. Disinfectant products composed either of didecyldimethylammonium chloride, quaternary ammonium carbonate or neutral enzymatic product. (e.g. Neodisher® mediclean) are also allowable.

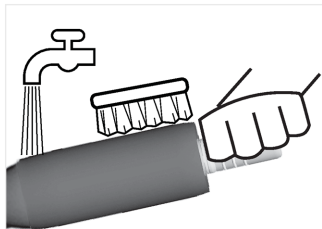


FIG. 7

5.2 Cleaning and disinfection

CAUTION

- Do not submerge in physiological liquid (NaCl) nor use saline solution to keep the device moist until it can be cleaned.
- Do not submerge in a cleaning bath.
- Do not clean in an ultrasonic cleaner.
- Do not spray any cleaning solution into the motor.
- Always ensure that the motor contacts are kept clean.

The external surface of the motor must be cleaned to remove impurities as follows [FIG. 7](#):

1. Hold the motor by the nose under tap water at 15°C - 38°C (59°F - 100°F) provided that the local tap water has a pH within the range of 6.5 - 8.5 and a chloride content below 100 mg/L. If the local tap water does not meet these requirements, use demineralized (deionized) water instead.
2. With the aid of a smooth flexible brush, clean the external surface of the motor using recommended cleaning products (e.g alkaline detergent or detergent-disinfectant pH 8-11).
3. Do not allow water to seep in the motor either by the nose or by the hose connector.
4. Dry the external surface of the motor with low linting textile impregnated with Spraynet®.

Nota: Automatic cleaning-disinfection can replace the above steps 2. If the automatic cleaner/disinfector does not provide efficient drying, follow the above step 4.

Washer-disinfector:

Carry out automatic cleaning- disinfection using an approved washer-disinfector which complies with ISO standard 15883-1.

Detergent and washing cycle:

Use an alkaline or detergent recommended for cleaning in a washer-disinfector for dental or surgical instruments (pH 8-11).

Recommended specifications for the thermo-disinfection cycle:

Phase	Parameters
Pre-cleaning	<45°C (113°F); ≥ 2 minutes
Cleaning	55°C — 65°C (131°F — 149°F); ≥ 5 minutes
Neutralization	≥ 2 minutes
Rinsing	Tap water, ≤30°C (86°F), ≥ 2 minutes cold water
Thermal Disinfection	Deminerlized water, 90°C — 95°C (194°F — 203°F), 5-10 minutes
Drying	18 — 22 minutes

CAUTION

Never rinse the devices to cool them.

CAUTION

If an automatic washer is used at the place of the washer/thermo-disinfector, respect the previous program for the Pre-cleaning, Cleaning, Neutralization and Rinsing phases. If the local tap water has a pH outside the range of 6.5-8.5 or if it contains more than 100 mg/l chloride (Cl-ion), do not dry the device inside the automatic washer but dry it manually with low linting textiles.

CAUTION

If automatic cleaner/disinfector does not provide efficient drying and/or if traces of humidity remain after drying, dry the external surface of the motor with low linting textile impregnated with Spraynet®.

5.3 Sterilization

CAUTION

- The quality of the sterilization is highly dependent on how clean the device is. Only perfectly clean devices may be sterilized.
- To improve the effectiveness of the sterilization, make sure the device is completely dry.
- Do not use a sterilization procedure other than the one described below.
- Only use dynamic air removal cycles: pre-vacuum or steam flush pressure pulse (SFPP) cycles.
- If the sterilization is required by national directives, use only dynamic sterilizers: do not use a steam sterilizer with a gravity displacement system. As with all instruments, following each sterilization cycle, including drying, remove the device to avoid an excessive exposure to heat which can result in corrosion.

5.3.1 Procedure

1. Pack the device in a packaging approved for steam sterilization.
2. Sterilize using steam, following dynamic air removal cycle (ANSI/AAMI ST79, Section 2.19), i.e. air removal via forced evacuation (ISO 17665-1, ISO/TS 17665-2) at 135°C (275°F), for 3 minutes or at 132°C (269.6°F) for 4 minutes. In jurisdictions where sterilization for prions is required, sterilize at 135°C (275°F) for 18 minutes.

The device sustains more than 1000 sterilisations.

The recommended parameters for the sterilization cycle are:

- The maximum temperature in the autoclave chamber does not exceed 137°C (278.6°F), i.e. the nominal temperature of the autoclave is set at 134°C (273.2°F), 135°C (275°F) or 135.5°C (275°F) taking into account the uncertainty of the sterilizer with regard to temperature.
- The maximum duration of the interval at the maximum temperature of 137°C (278.6 °F) is in accordance with national requirements for moist heat sterilization and does not exceed 30 minutes.
- The absolute pressure in the chamber of the sterilizer is comprised in the interval 0.07 bar to 3.17 bar (1 psia to 46 psia).
- The rate of change of temperature does not exceed 15°C/min (59°F/min) for increasing temperature and -35°C/min (-31°F/min) for decreasing temperature.
- The rate of change of pressure does not exceed 0.45 bar/min (6.6 psia/min) for increasing pressure and -1.7 bar/min (-25 psia/min) for decreasing pressure.
- No chemical or physical reagents are added to the water steam.

5.4 Lubrification

CAUTION

Do not spray any lubricant or cleaning-lubricant solution into the motor [FIG. 8](#).

5.5 Packing and storage

Storage conditions



Temperature range:

0°C — +40°C (+32°F — +104°F)



Relative humidity range:

10% — 80%



Air pressure range:

650 hPa — 1060 hPa



Keep away from rain

The device must be stored inside the sterilization pouch in a dry and dust free environment. The temperature must not exceed 55°C (131°F). If the device will not be used for 7 days or more after the sterilization, extract the device from the sterilization pouch and store it in the original package. If the device is not stored in a sterilization pouch or if the pouch is no longer sterile, clean, dry and sterilize the device before using it.

CAUTION

If the medical device has been stored refrigerated, allow it to warm up to room temperature prior to its use.

CAUTION

Comply with the expiration date of the sterilization pouch which depends on the storage conditions and type of packaging.

5.6 Servicing

Bien-Air Dental SA recommends that the user has his or her dynamic devices checked or serviced after 5000 processing cycles or 5 years.





CAUTION

Never disassemble the device. For all modification and repair, contact your regular supplier or Bien-Air Dental service centre.

6 Transport & disposal

6.1 Transport

Transport conditions

	Temperature range:	-20°C — +50°C (-4°F — +122°F)
	Relative humidity range:	5% — 80%
	Air pressure range:	650 hPa — 1060 hPa
	Keep away from rain	

6.2 Disposal



The disposal and/or recycling of materials must be performed in accordance with the legislation in force.



This device must be recycled. Electrical and electronic equipment may contain dangerous substances which constitute health and environmental hazards. The user must return the device to its dealer or establish direct contact with an approved body for treatment and recovery of this type of equipment (European Directive 2012/19/ EU).

7 General information

7.1 Terms of guarantee

Bien-Air Dental SA grants the operator a guarantee covering all functional defects, material or production faults.

The warranty period is:

- 36 months from the date of invoicing.

In the event of justified claim, Bien-Air Dental or its authorised representative will fulfil the company's obligations under this guarantee by repairing or replacing the product free of charge.

Any other claims of any kind whatsoever, particularly claims for damage or injury and the consequences thereof resulting from:

- Excessive wear and tear
- Infrequent or improper use
- Failure to observe the servicing, assembly or maintenance instructions
- Damage caused by unusual chemical, electrical or electrolytic influences
- Faulty air, water or electrical connections



CAUTION

The warranty becomes null and void if damage and its consequences result from incorrect servicing or modification by third parties not authorized by Bien-Air Dental SA. Warranty requests will only be taken into consideration if the product is accompanied by a copy of the invoice or delivery note. The following information must be clearly indicated: purchase date, product reference and serial number.

7.2 Cable compatibility

REF	Micromotor	REF Cable
1600755-001	MOT MX-i LED	1601069-001
1600825-001	MOT MX-i	1601881-001
1601008-001	MOT MX-i LED 3RD GEN	1601009-001

7.3 References

REF	Legend
1600825-001	MX-i micromotor
1601008-001	MX-I LED 3RD GEN micromotor
1600755-001	MX-i LED micromotor
1600036-006	Spraynet®, cleaning spray 500 ml, box of 6 cans
1601069-001	MX-i LED, cable
1601009-001	MX-I LED 3RD GEN cable
1600881-001	MX-i LED cable (3 meters length)
1300967-010	O-RING 8.1x0.73



Bien-Air Dental SA

Länggasse 60 Case postale 2500 Bienne 6 Switzerland

Tel. +41 (0)32 344 64 64 Fax +41 (0)32 344 64 91

dental@bienair.com

Other addresses available at

www.bienair.com

EC REP Bien-Air Europe Sàrl

19-21 rue du 8 mai 1945

94110 Arcueil

France